

ALPHKUT® BRAVO [AC-400CMB]

- The PCS System gives an audio-visual alarm & stop the RF output for safe procedures
- Power, Modes Settings are accessible through touch keyboard
- The System memory offers 99 user-defined programs, plus a preset SUPPLIED program
- Peak power system for optional support during the initial cutting in under water surgery
- Newly develop bipolar cutting mode | more safety, current flow only at the target tissue.
- Extremely easy to operate, thanks to the supported touch keyboard
- Monopolar cutting & coagulation done by using both hand & foot switch
- Isolated Rf output to the generator significantly reduces the risk of alternate site burn.
- Automatic output adjustment dependent on tissue impedance
- Latest advanced micro controller based technology used which is upgradable.
- Memorization of the last operating condition.
- Lesser charring and drag
- Bipolar cut function for sealing & cutting simultaneously.
- Permanent & safe sealing up to 7mm.
- Minimal lateral thermal spread.
- Resection can be perform with saline also.



DEVICE

Part No.	Description
AC-400CMB	Alphkut® Bravo [Complete Monopolar & Bipolar] [Electrosurgical Unit] With All Its Standard Accessories 400W (Monopolar, Bipolar Bi-cut, Bi-tur & Sealer)

SUPPLIED ACCESSORIES

Part no.	Description
SSP-02	SS Patient plate
PPC-02	Patient plate cable
FSC-03T-CMB	Tri-Foot switch with toggle
MCH-01	Monopolar chuck handle
MC-03	Mains Cord
NLPG-05-30-31-SA	Lap Bi-Clamp Grasper
NLCP-05-30-31-SA	Lap Bi-Clamp Cutter
BCS-04-03-VMF	Bipolar cable

TECHNICAL DATA

Power supply	170 - 260 VAC
Maximum power	400W
Classification	Protection Class 1, type Bf
Dimension (w x h x d)	325 x 125 x 360mm
Weight	8kg

TECHNICAL DATA

Monopolar cut

Pure cut	400 Watt	300 OHM	1.5 CF
Blend I	250 Watt	300 OHM	2.0 CF
Endocut	100 Watt	300 OHM	1.5 CF

Monopolar coag

Spray	120 Watt	500 OHM	7.0 CF
Fulgurate	150 Watt	500 OHM	6.0 CF
Soft	120 Watt	500 OHM	5.0 CF

Bipolar

Micro	70 Watt	100 OHM	1.5 CF
Macro	90 Watt	100 OHM	1.5 CF
Bi-Seal	4 Level		
P: pulsating power	P1,P2,P3,P4	100 OHM	1.5 CF
C: continuous power	C1,C2,C3,C4	100 OHM	1.5 CF

Bi-Cut

Cut (10-90)	200 Watt	100 OHM	1.5 CF
-------------	----------	---------	--------


Bi-tur Cut

Pure Cut	300 Watt	100 OHM	1.5 CF
Blend cut	250Watt	100 OHM	1.7 CF

Bitur Coag

Spray	120Watt	100 OHM	4.5 CF
Fulguration	150Watt	100 OHM	4.0 CF
Dessication	120Watt	100 OHM	3.5 CF


SS Patient plate
SPP-02

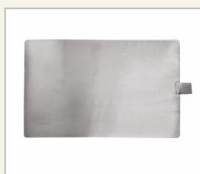

Patient plate cable
PPC-02



Tri-foot switch with toggle
FSC-03T-CMB

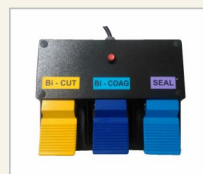

Mains cord
MC-03



Monopolar chuck handle
MCH-01



Lap Bi-Clamp Grasper
NLPG-05-30-31-SA



Lap Bi-Clamp Cutter
NLCP-05-30-31-SA



Bipolar cable
BCS-04-03-VMF



*Open bipolar tweezer
OBT-07



*Silicon Patient plate
SPP-02



*Open Bipolar Clamp
OBPR-08



*Open Bipolar Clamp
BCSO-04-03-3MF



Further Accessories available on request