


**ALPHSORCE® LED LIGHT [LSL-100]**

- Long life > 20000 hrs and more
- Good Lumen output ( Superior to 250W of halogen)
- Instant on/off (no warm up or cool down)
- No change in color temperature over entire range of colour spectrum.
- Does not emit ultraviolet and infrared light. Reduced heat transmission to the Endoscope
- LCD Display : Brightness Control
- Consumes very less power
- No annual bulb change is necessary

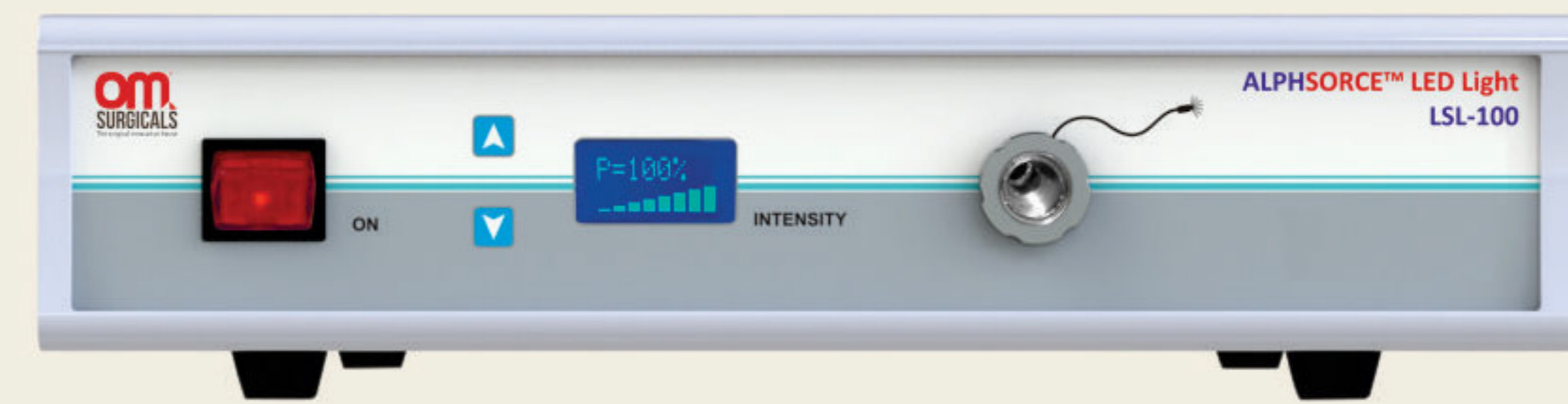

**DEVICE**

Part No.	Description
LSL-100	Alphsorce™LED Light L100 for Laparoscopy with Light cable socket,Storz

**SUPPLIED ACCESSORIES**

Part no.	Description
MC-03	Mains cord

**TECHNICAL DATA**

Lamp	LED LSL-100
Color Temperature	6000°k
Light guide connection	Karl Storz type, others on request
Lumen Maintenance	90% After 10000hrs
Lumen output	750 lumens
Ambient temperature	+10°C to +40°C in operating mode
Cooling Method	Forced Air cooling
Power Supply	100-240 AC
Frequency	50 Hz
Maximum Power	25W
Classification	Protection Class I, Type BF
Dimensions (W x H x D)	300 x 30 x 255 mm
Weight	5.5 Kgs.

Mains cord  
MC-03



\*Fibre optic cable Armoured  
FOC-523-SS



\*Fibre optic cable Non-armoured  
FXN-523-SS

